




User Manual

2BCURV BELLY & LUMBAR MASSAGER

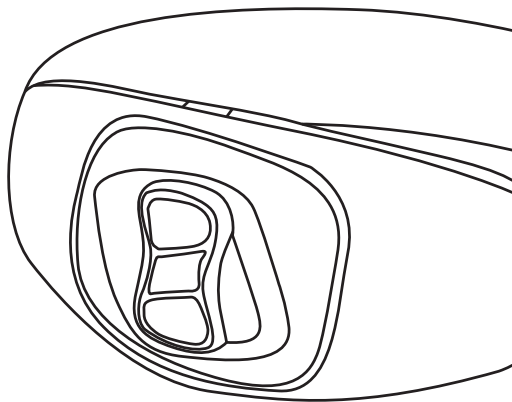

MN-VM02.

Please read and keep this manual carefully before using the product.

CONTENT

1. Product Components	02
2. How to Use 2BCurv	03
3. How to Use 2BCurv - Remote Control	04
4. Product Specifications	05
5. Safety Instructions	06
6. Troubleshooting	08

PRODUCT COMPONENTS

Figure 1


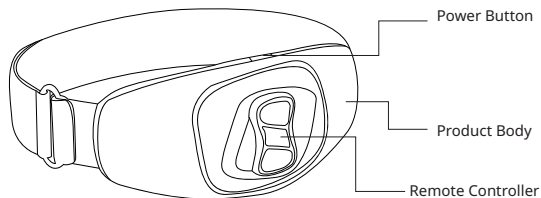



Figure 2



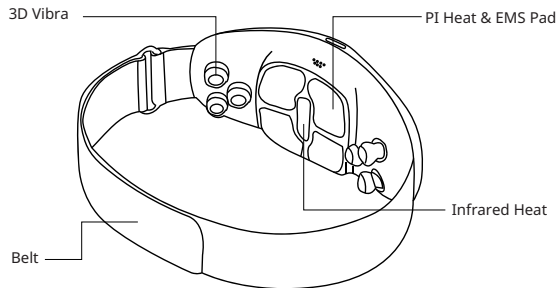

HOW TO USE 2BCURV


1. Put the belt around the area to be massaged, tighten the belt. Making sure the belt is fixed on the body comfortably.
2. Press and hold the power button to start 2BCurv. Follow the remote-control operation in next page to operate the device.
3. After completion of use, you are recommended to drink a glass of water.
4. After used, it is recommended to store it tidily.
5. Please wear cotton clothes while operating 2BCurv.

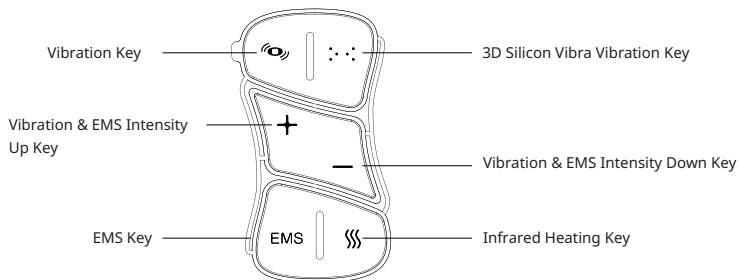CONTROL PANEL



HOW TO USE 2BCURV - REMOTE CONTROL

1. Press and hold the Power button to start the device. The Power button will light up blue when in operation. Press and hold the Power button again to turn off the device. When charging, the Power button will light up in pulsing blue. When device is fully charged, the light will stay constant blue.
2. **Vibration Key**  : Pressing Vibration key on remote activates vibration mode 1 by default. At the same time the remote key will light up white. There are 5 modes available. Short press this key again to circle through mode 2, mode 3, mode 4 and mode 5 then back to mode 1 and so on.
3. **3D Silicon Vibra Vibration Key**  : Press this key to activate the 3D Silicon Vibra function mode 1 by default. At the same time the remote key will light up white. There are 4 modes available. Short press this key again to circle through mode 2, mode 3 and mode 4 then back to mode 1 and so on.
4. **+ - Key**: Control the intensity of vibration and EMS, 4 intensities available for vibration mode. 16 intensities available for EMS mode. Short press the key for adjustment. Default is intensity level 1. Button will light up white when pressed.
5. **EMS Key**: Press this key to activate the EMS function mode 1 by default. At the same time the remote key will light up white. There are 6 modes available. Short press this key again to circle through mode 2, mode 3, mode 4, mode 5 and off.

6. **Infrared Heating Key**  : Press this key to activate the Infrared heating function mode 1 by default. At the same time the remote key will light up white. There are 3 modes available. Short press this key again to circle through mode 2 and off.

Remark: Vibration and EMS + Heating mode is mutually exclusive. They cannot work together due to battery limitation.

PRODUCT SPECIFICATIONS

Product Name	2BCurv Belly and Lumbar Massager
Model No.	MN-VM02
Rated Power	10W
Battery	1500mAh x 3pcs
Power Input	5V DC, 2A
Timer	15 Mins
Charging Type	Type C

1. Unplug cable from this product after use and before cleaning and maintenance.
2. After each usage, use soft cloth to clean and cover the massager before storing.

3. After cleaning, use dry towel to wipe.
4. Keep the product on a proper dry storage area and keep away from water and wet places. Do not use tools to disassemble the massager for cleaning.
5. Do not put under direct sunlight or rain; put in cool and dry place after use.

SAFETY INSTRUCTIONS

PROHIBITION

1. Prohibited from using this product while bathing.
2. Avoid direct contact with water.
3. Do not use in temperature less than 0 ° C
4. Do not use while in hunger and or after you have just completed eating.
5. Do not massage on injured skin.
6. Prohibited to pull out the electric conductor.

7. In case of any trouble or abnormal phenomenon/operation, stop using immediately. Contact MYNT customer service to troubleshoot.
8. Do not use in wet environment to avoid water from penetrating into the massager that will cause danger.

WARNINGS

1. If you are in medical treatment or having any cases as below, please be careful in using this product. Get advice from your doctor before using: Diabetes, Stroke, Liver diseases, Tumor, Skin diseases, Infectious diseases, Adenitis, Cardiovascular diseases (Hypertension.)
2. If you have an open injury or illness, please consult your doctor before use.
3. Not suitable for pregnant women and children, women in menstruation should consult doctor before use.
4. To avoid any unforeseen circumstances, people with nervous disorder are advice not to operate the massager alone.
5. Children and people with mental deficiency must be monitored and supported by adult in use.
6. When the massager is in operation, do not let children play besides.

7. Do not let children under 12 years old to use.
8. Heat generator exists inside the product, massager surface may be hot. For people not sensitive to heat, please be very careful when using this product for prolong session.

TROUBLESHOOTING

FAILURE	POSSIBLE CAUSES	SOLUTIONS
Product cannot operate	Low battery	Charge the device
Why all function cannot be used together	Shaking + 3D one group; EMS + Heat + 3D another group	Design based on product features and battery limitation
No power	Power button is not switched on	Press power button
Not vibrating	Vibrating function is not started	Start vibrating function

Vibration not working	Heat or EMS is working	Turn off EMS or Heat function and swap to Vibration mode
No Heating	Vibration is working	Turn off Vibration function and swap to Heat function
No EMS	Vibration is working	Turn off Vibration function and swap to EMS function
Unable to charge	USB charging cable is not inserted properly	Insert USB charging cable correctly

Contact MYNT Customer Service if problem still persist.



SCAN ME

REGISTER E-WARRANTY in 4 SIMPLE STEPS

Please register your warranty online



01

Scan QR Code or
log in to
<https://www.mynt.com.my/e-warranty.html>



03

Register Product
Warranty



02

Attach Product
Invoice



04

Receive Confirmation
Email