

SKINRUB™ LITE OPERATING MANUAL

PRODUCT PARAMETERS

Product name : Skinrub Lite Massager
Product model : ZD-JS0201
Implementation standard : Q/HNZD 013-2020
Rated voltage : 5V =
Rated power : 5W

REMOVE BATTERY

When removing the battery, be sure to protect your eyes, hands and any parts that may cause personal injury.

Please follow these steps:

1. Use a screwdriver to open the casing and take out the internal components.
2. Use a soldering iron to melt the soldering point, disconnect, and pry out the battery (non-professionals should not dismantle it).

zdeer

mynt



SCAN ME

REGISTER E-WARRANTY in 4 SIMPLE STEPS

Please register your warranty online



01

Scan QR Code or
log in to
<https://www.mynt.com.my/e-warranty.html>



03

Register Product
Warranty



02

Attach Product
Invoice




04

Receive Confirmation
Email

SKINRUB™ LITE

PRODUCT DESCRIPTION

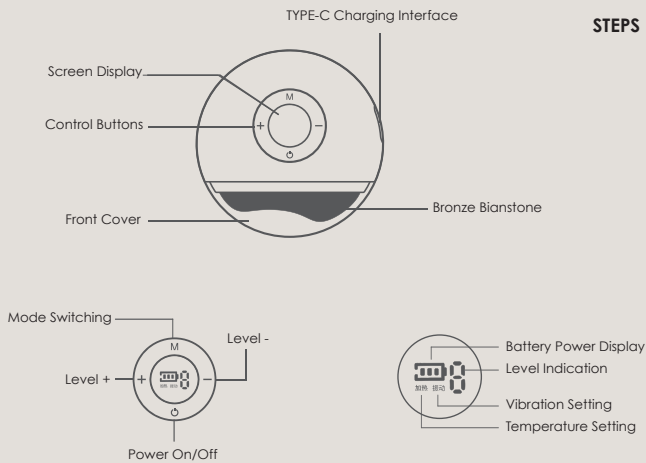

STEPS FOR USAGE

Power on

"Press and hold the power button for 3 seconds, the buzzer will sound, the display will light up, and the product will enter the power-on state. (The default setting is the heating interface when turned on)

Temperature setting

In the heating setting interface, short press the "+" key. The buzzer will sound and the level will increase by one. There are 9 heating levels in total (level 0 is standby mode). After reaching level 9, short pressing will not increase the level further. Short press the "-" key, buzzer will sound, and the level will be reduced by one level. After reaching level 0, short pressing will not reduce the level

Vibration setting

In the vibration setting interface, short press the "+" key. The buzzer will sound, and the level will increase by one. There are 9 vibration levels in total (level 0 is standby mode). After reaching level 9, short pressing will not increase the level further. Short press the "-" key, the buzzer will sound once, and the level will be reduced by one level. After reaching level 0, short pressing will not reduce the level (temperature and vibration functions can be used independently)

Shut down

Press and hold the power button for 3 seconds. The buzzer will sound, the display will go out, and the product will enter the shutdown state. (If there is no operation within 15 minutes, the product will automatically shut down) vibration functions can be used independently)

Note: The product cannot be turned on and used while charging. It must be turned on by unplugging the power source

Warning: The appliance has a heating surface. People who are not sensitive to heat should pay attention when using it

PRODUCT USE SCHEMATIC DIAGRAM

The recommended order of use is head, neck, back (thoracic, lumbar, sacral), chest, abdomen, upper limbs, and lower limbs.




TOXIC AND HAZARDOUS SUBSTANCES CONTENT TABLE

Part Name	Toxic and harmful substances or elements					
	Lead (Pb)	Hg (Hg)	Cadmium (Cd)	Hexavalent Chromium (Cr(VI))	Polybrominated biphenyls (PBB)	Polybrominated diphenyl ethers (PBDE)
Shell	○	○	○	○	○	○
Circuit Board	×	○	○	○	○	○
Power Cable	○	○	○	○	○	○
<p>10 This form is prepared in accordance with the provisions of SJ/T11364</p> <p>O: Indicates that the content of harmful substances in the material of this part meets the limit requirements specified in GB/T26572</p> <p>X: Indicates that the content of hazardous substances in the material of this part exceeds the requirements of GB/T26572. Certificate</p> <p>Note: Parts containing toxic and hazardous substances or elements cannot be replaced by toxic and hazardous substances or elements due to limitations of global development levels</p>						

CONTRAINDICATIONS

- Critical illness, such as severe heart disease, renal failure, cirrhosis, ascites, severe general edema and etc. are prohibited from use.
- People with highly sensitive skin, severe allergies, and patients with contagious skin diseases are prohibited from using this product.
- Do not use this product on areas with boils, ulcers, sores, rashes, or unexplained masses on the body surface.
- People with bleeding tendencies, such as advanced diabetes, severe anemia, leukemia, aplastic anemia and thrombocytopenia and etc., are prohibited from using this product.
- People with low blood pressure, hypoglycemia, excessive weakness, nervousness, or a fear of pain are prohibited from using this product.
- This product is prohibited for use by the elderly, those who are extremely weak or emaciated, or those who have lost skin elasticity due to long-term illness.
- Do not use on painful areas of acute sprains, trauma, or fractures.
- Do not use this product if you are overly thirsty, hungry, drunk, tired, or stay up too late.
- Do not use on the lower abdomen during menstruation. Pregnant or breastfeeding women are not allowed to use it.

PRECAUTIONS AND SAFETY

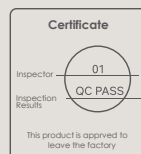

- Pay attention to the temperature of the copper Bianstone when using it, and choose a suitable temperature level based on what your skin can tolerate. Do not overheat to prevent burns.
- If there is any odor residue after using this product with essential oils, wipe it with a clean damp cloth to keep the product clean. Disinfect with 75% alcohol before using it on others to avoid infection.
- Do not move the instrument by carrying the power cord.
- This product cannot be dropped or hit from a high altitude. Handle it with care, as it is fragile.
- This machine is not waterproof. Do not immerse it in water or splash it with water.
- It is strictly prohibited to use harsh chemicals to clean the machine body, power supply and button cover.
- The appliance has a heating surface, so people who are not sensitive to heat must pay attention when using it.

Note: This product cannot be used as a substitute for medicines or medical devices

OPERATING MANUAL

Product model ZD-JS0201 Rated Voltage 5.0V= Rated Power 5W Implementation Standard Q/HNZD 013-2020 Brand Operation Zuodian Industrial (Hubei) Co., Ltd. Operation Address 1st floor, Block C, Jinglun Park, No. 70 Guanggu Avenue, Wuhan City, Hubei After-Sales Service Unit Zuodian Industrial (Hubei) Co., Ltd. After-Sales Service Phone Number 400-035-6661 Manufacturer Henan Zuodian Health Industry Co., Ltd. Production address Zuodian Science and Technology Park, No. 1 Jiankang Road, Jiankang Road, Changshu Road, Changge, Xuchang, Henan, China. This product is not a substitute for medical devices or medicines.